

*EURL for Cereals and Feeding stuff  
National Food Institute  
Technical University of Denmark*

## **Validation Report 16**

**Determination of pesticide residues in wheat, barley and rice  
by GC-MS/MS and LC-MS/MS**

**(QuEChERS method)**

**Susan Strange Herrmann**

**Gitte Andersen**

**Mette Erecius Poulsen**

**December 2014**

**CONTENT:**

<i>1. Introduction</i> .....	3
<i>2. Principle of analysis</i> .....	3
<i>3. Validation design</i> .....	3
<i>4. Chromatograms and calibration curves</i> .....	4
<i>5. Validation parameters</i> .....	8
<i>6. Criteria for the acceptance of validation results</i> .....	9
<i>7. Results and discussion</i> .....	10
<i>8. Conclusions</i> .....	11
<i>9. References</i> .....	11
<i>Appendix 1a. MRM transitions GC-MS/MS</i> .....	12
<i>Appendix 1b. MRM transitions for LC-MS/MS</i> .....	13
<i>Appendix 2. Recoveries, repeatability (RSD<sub>r</sub>), internal reproducibility (RSD<sub>r</sub>) and Limit of Quantification (LOQ) for pesticides validated on three cereal commodities, wheat, barley and rice.</i> .....	14
<i>Appendix 3. Recoveries, repeatability (RSD<sub>r</sub>) and Limit of Quantification (LOQs) for pesticides validated on wheat.</i> .....	16
<i>Appendix 4. Recoveries, repeatability (RSD<sub>r</sub>) and Limit of Quantification (LOQs) for pesticides validated on barley</i> .....	18
<i>Appendix 5. Recoveries, repeatability (RSD<sub>r</sub>) and Limit of Quantification (LOQs) for pesticides validated on rice</i> .....	20
<i>Appendix 6: Principles of the QuEChERS method for cereal extraction</i> .....	22

## 1. Introduction

This report describes the validation of the QuEChERS method combined with GC-MS/MS and LC-MS/MS. The method was sought validated for 34 pesticides in wheat, rice and barley. The QuEChERS method is an extraction method which has been developed to be Quick, Easy, Cheap, Efficient, Rugged and Safe. The method is most commonly used on fruit, vegetables and cereals<sup>1</sup>.

## 2. Principle of analysis

**Sample preparation:** The samples is milled with a sieve at 1 mm.

**Extraction:** The sample is shaken and a salt and buffer mixture is added and the sample is shaken again.

**Clean-up:** After centrifugation the supernatant is transferred to a clean tube and put in -80 degree freezer. When the extract is almost thawed it is centrifuged and the supernatant is transferred to a tube containing PSA and MgSO<sub>4</sub>. After shaking and an additional centrifugation step the final extract is diluted 1:1 with acetonitrile to obtain the same matrix concentration as in the matrix matched calibration standards.

**Quantification and qualification:** The final extract is analysed by GC/MS/MS and LC-MS/MS.

**GC-MS/MS:** The pesticide residues are separated on a DB5-MS column and analysed by triple quadrupole operating in the multiple reaction monitoring mode (MRM) with electron energy at 70 eV, source temperature at 180°C and transfer line at 250°C. The injection volume was 4 µl. For each pesticide two sets of precursor and product ions were determined. One for quantification and one for qualification. The MRM transitions for the pesticides and degradation products are given in **Appendix 1a**.

**LC-MS/MS:** The pesticide residues are separated on a reversed-phase column and detected by tandem mass spectrometry (MS/MS) by electrospray (ESI). The validation includes pesticides determined with both positive and negative ESI. <sup>13</sup>C<sub>6</sub>-carbaryl was used as internal standard but was not used for the quantification. All pesticides were detected in the MRM mode. For each pesticide precursor ion and 2 product ions were determined. One product ion for quantification and one for qualification. The MRM transitions for the pesticides and degradation products sought validated are given in **Appendix 1b**.

## 3. Validation design

The method was south validated for 34 pesticides or degradation products in wheat, rice and barley, see **Table 1**. The validation was performed on 5-6 replicates on each cereals commodity at each of


the three spiking levels; 0.01, 0.02 and 0.1 mg/kg. A blank sample of each cereal commodity was included.

**Table 1. Pesticides included in the recovery experiments.**


<b>Pesticides included in recovery experiments</b>		
2-hydroxypropoxycarbazone	Dichlobenil	Isocarbophos
4-trifluoromethylnicotinic acid (TNFA)	Dimethenamid	Meptyldinocap
Acibenzolar-S-Methyl	Ethalfuralin	Metobromuron
Anthraquinone	Ethofumesate	Nicosulfuron
Asulam	Ethoxysulfuron	Propachlor
AvermectinB1a (Abamectin)	Etoxazole	Propoxycarbazone-sodium salt
Beflubutamid	Fenpyroximate	Pyraflufen-Ethyl
Bifenazate	Flumioxazin	Rotenone
Biphenyl	Foramsulfuron	Spirotetramat
Butralin	Forchlorfenuron	Tetramethrin
Carfentrazone-Ethyl	Imazosulfuron	Trinexapac-ethyl
Chlorotoluron		

#### 4. Chromatograms and calibration curves


The calibration curve is determined by the analysis of each of the analysts at least 4 calibration levels, i.e. 0.003, 0.01, 0.033 and 0.1 µg/ml. The calibration curves were in general best fitted to a linear curve. The quantification was performed from the mean of two bracketing calibration curves. The majority of the correlation coefficients (R) were higher or equal to 0.99. Examples of chromatograms obtained when analysing the extracts by GC-MS/MS are presented in **Figure 1-4**. Examples of calibration curves for LC-MS/MS are presented in **Figure 5-8**.

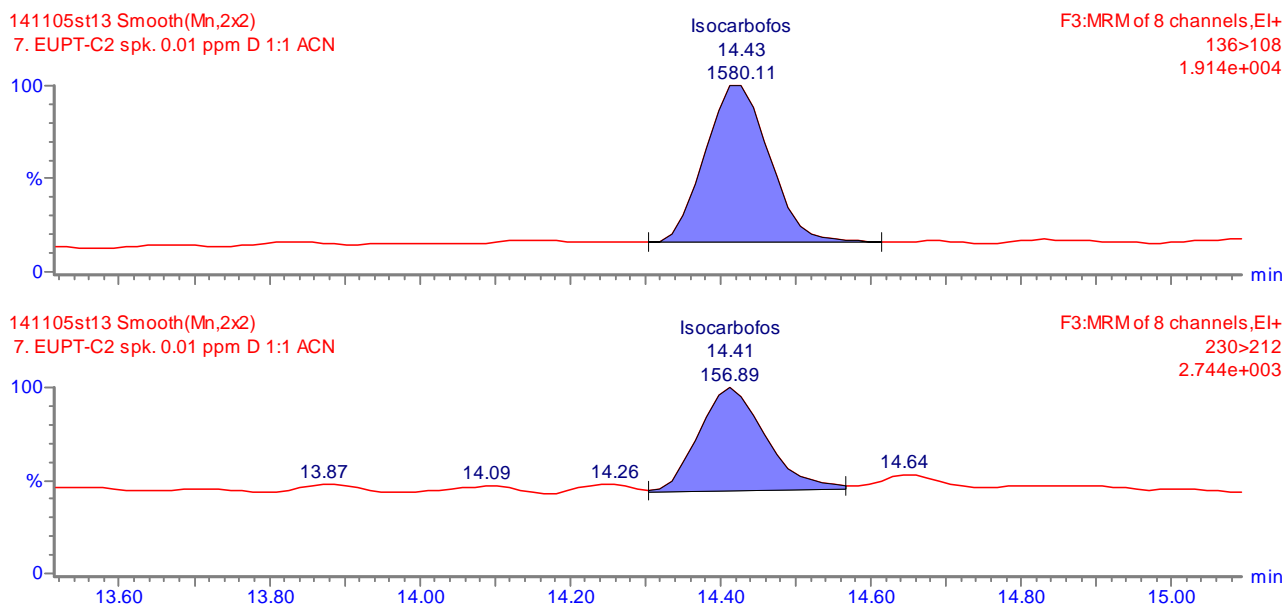


**Figure 1:** Examples of GC-MS/MS chromatograms for isocarbofos in wheat obtained when analysing extract spiked with 0.01 mg/kg (two MRM transitions are shown for each pesticide).

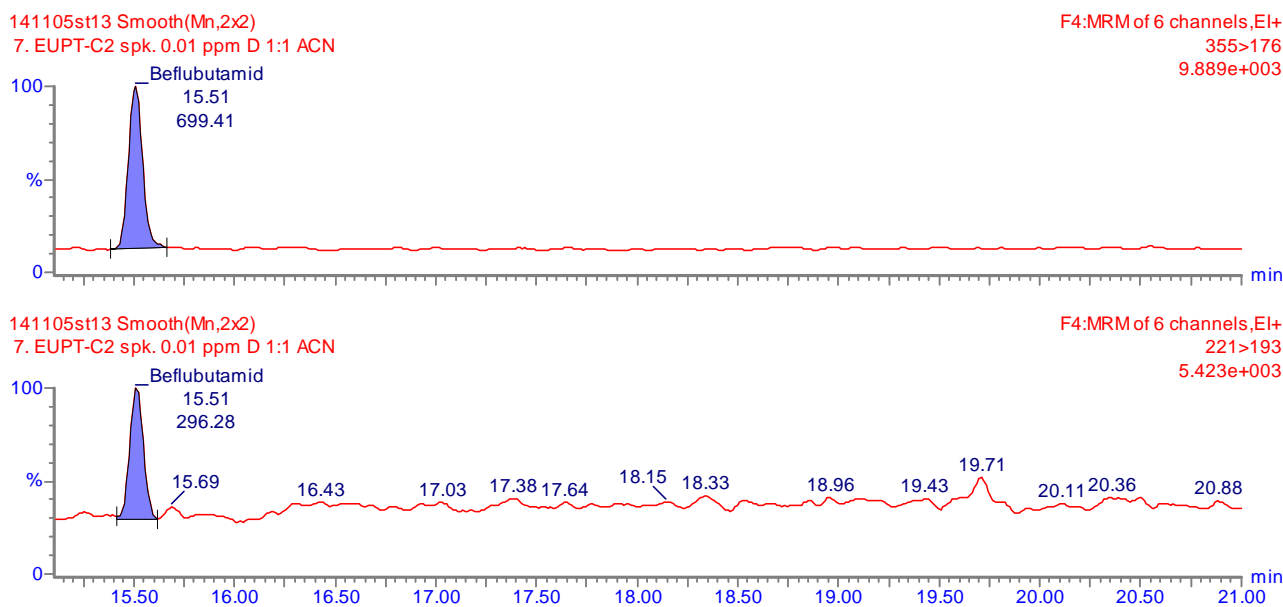

**Figure 2:** Examples of GC-MS/MS chromatograms beflubutamid obtained when analysing extract spiked with 0.01 mg/kg.

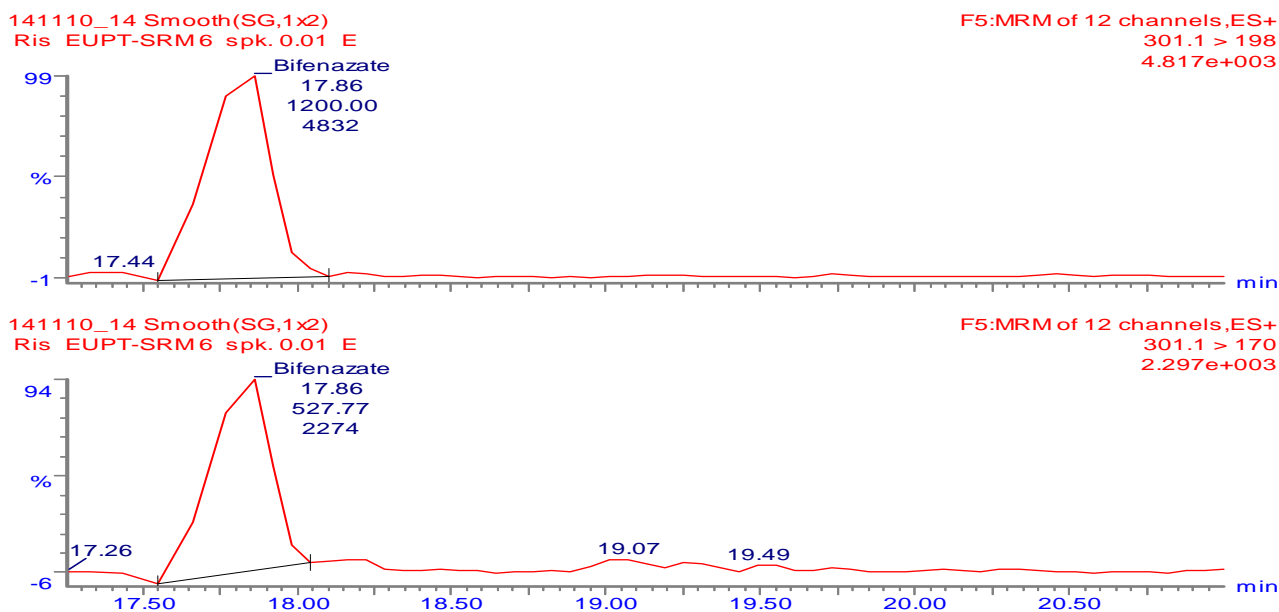


**Figure 3:** Examples of LC-MS/MS chromatograms bifenazate in rice obtained in positive mode when analysing extract spiked with 0.01 mg/kg (two MRM transitions are shown for each pesticide).

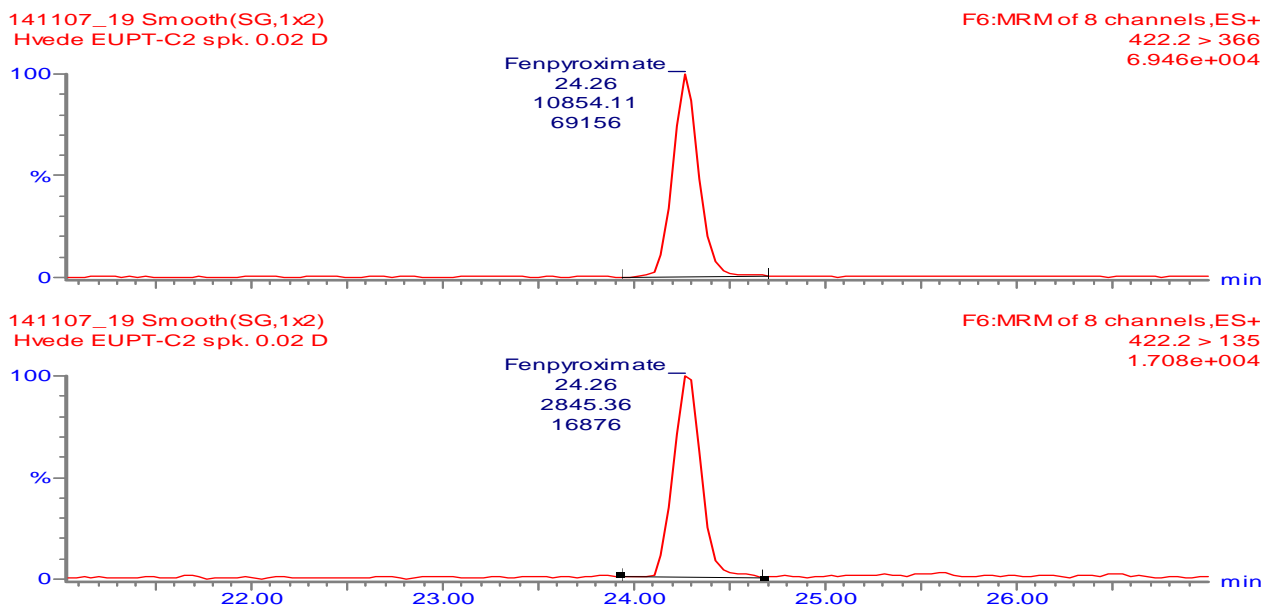


**Figure 4:** Examples of LC-MS/MS chromatograms fenpyroximate in rice obtained when analysing extract in positive mode spiked with 0.02 mg/kg.

Compound name: Isocarbofos  
 Correlation coefficient:  $r = 0.992100$ ,  $r^2 = 0.984262$   
 Calibration curve:  $656094 * x + 26.4734$   
 Response type: External Std, Area  
 Curve type: Linear, Origin: Exclude, Weighting: 1/x, Axis trans: None

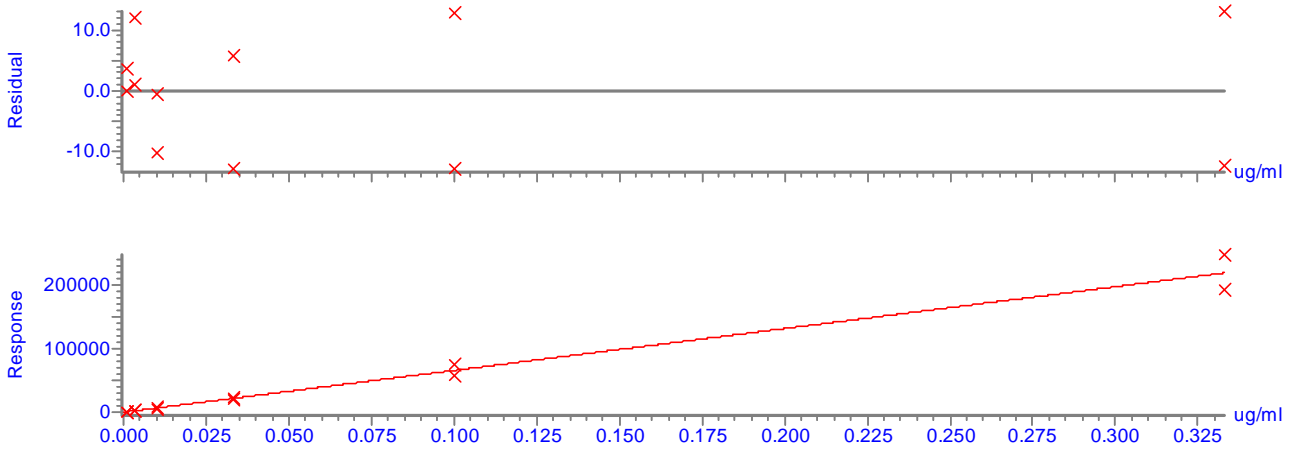


**Figure 5.** Examples of GC-MS/MS calibration curves for isocarbofos matrix matched with wheat (concentrations from 0.001-0.333  $\mu\text{g/ml}$ )

Compound name: Beflubutamid  
 Correlation coefficient:  $r = 0.996135$ ,  $r^2 = 0.992285$   
 Calibration curve:  $272220 * x + -4.79738$   
 Response type: External Std, Area  
 Curve type: Linear, Origin: Exclude, Weighting: 1/x, Axis trans: None

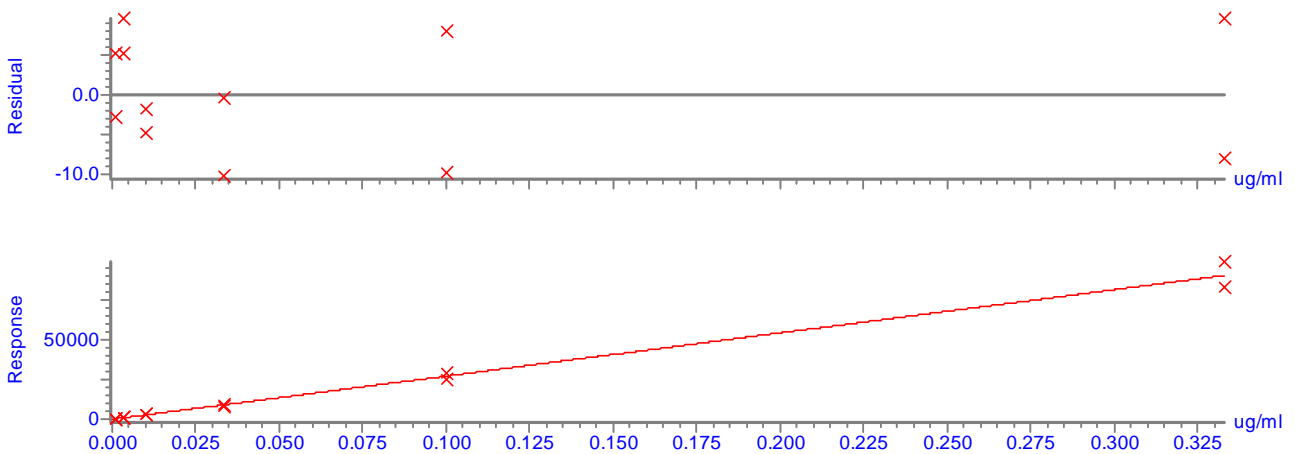


**Figure 6.** Examples of GC-MS/MS calibration curves for beflubutamid matrix matched with wheat (concentrations from 0.001-0.333  $\mu\text{g/ml}$ ).

Compound name: Bifenazate  
 Correlation coefficient:  $r = 0.999470$ ,  $r^2 = 0.998941$   
 Calibration curve:  $734791 * x + 71.3433$   
 Response type: External Std, Area  
 Curve type: Linear, Origin: Exclude, Weighting: 1/x, Axis trans: None

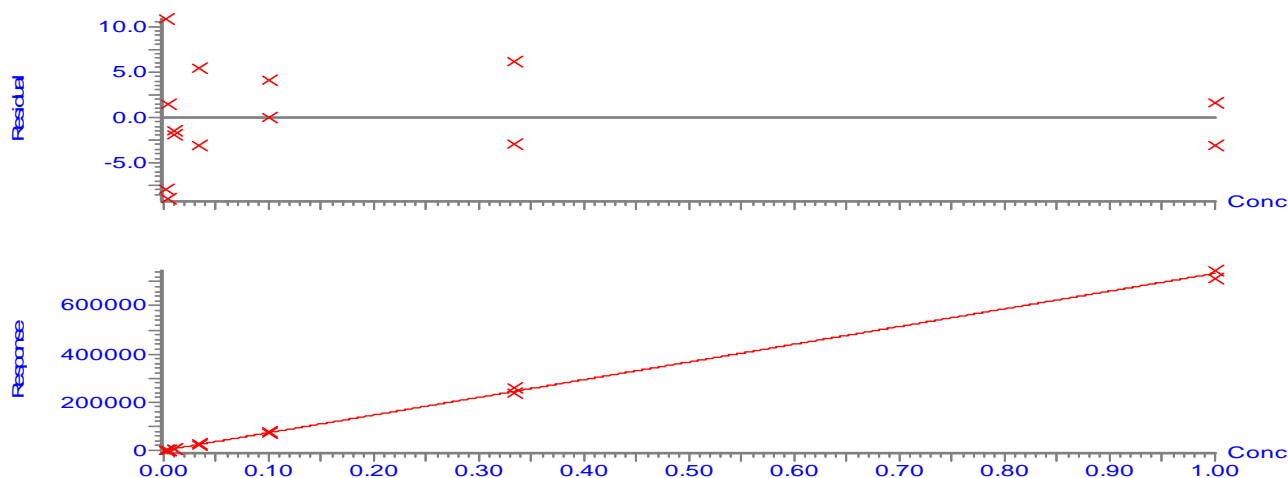

**Figure 7.** Examples of LC-MS/MS calibration curves for bifenazate matrix matched with wheat (concentrations from 0.001-1.0  $\mu\text{g/ml}$ ).

Compound name: Fenpyroximate  
 Correlation coefficient:  $r = 0.998700$ ,  $r^2 = 0.997401$   
 Calibration curve:  $2.54129\text{e}+006 * x + 253.549$   
 Response type: External Std, Area  
 Curve type: Linear, Origin: Exclude, Weighting: 1/x, Axis trans: None

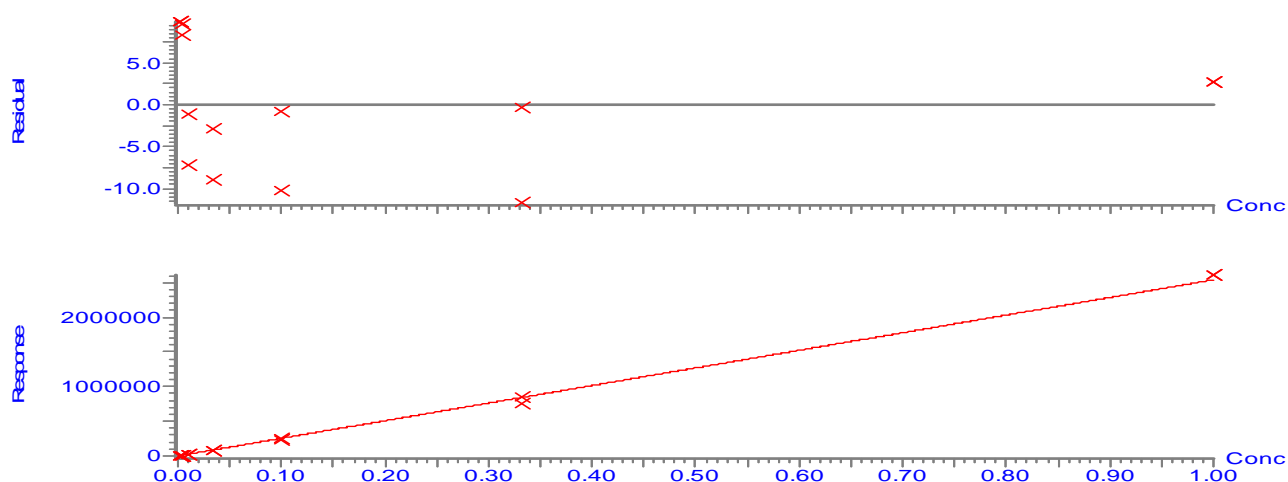

**Figure 8.** Examples of LC-MS/MS calibration curves for fenpyroximate matrix match with rice (concentrations from 0.001-1.0  $\mu\text{g/ml}$ ).

## 5. Validation parameters

### Precision – repeatability and internal reproducibility

Repeatability was calculated for all pesticides and degradation products on all three spiking levels (0.01 mg/kg, 0.02 mg/kg and 0.1 mg/kg), both for the individual cereal commodities and for the all commodities altogether. Repeatability is given as the relative standard deviation on the result from two or more analysis at the same sample, done by the same technician, on the same instrument and



within a short period of time. The internal reproducibility is calculated for the all the cereal commodities only, because the individual cereal type is analysed on one occasion only. Internal reproducibility is relative standard deviation on results obtained under reproducibility conditions, with the same method on the same sample by different operators within a larger period of time.

Repeatability ( $RSD_r$ ) and internal reproducibility ( $RSD_R$ ) in this validation was calculated from the 5-6 replicate determinations. Repeatability were calculated as given in ISO 5725-2<sup>2</sup>.

#### **Accuracy – Recovery**

The accuracy was determined from recovery studies in which samples were spiked at three concentration levels (0.01 mg/kg, 0.02 mg/kg and 0.1 mg/kg) with the relevant pesticides, isomers and degradation products.

#### **Robustness**

The QuEChERS method has, in connection with the development of the method, been shown to be robust by Anastassiades et al. 2003<sup>1</sup>.

#### **Limit of quantification, LOQ**

The quantification limits (LOQ) was determined as the lowest spike level for which the acceptance criteria (se Section 6) was meet.

The obtained results including recovery,  $RSD_r$ ,  $RSD_R$  and limit of quantification (LOQ) are presented in appendix 2 for the pooled results obtained for all three types of cereal and in appendix 3 to 5 for the individual cereal types; wheat, barley and rice

### **6. Criteria for the acceptance of validation results**

For the pesticides to be accepted as validated the following criteria for precision and trueness must to be fulfilled:

1. The relative standard deviation of the repeatability should be  $\leq 20\%$ <sup>3</sup>.
2. The average relative recovery must be between 70 and 120%<sup>3</sup>.

If the above mentioned criteria have been meet, the quantification limits, LOQs have been calculated.

## 7. Results and discussion

Initially acibenzolar acid, chloropicrin and dithianon were to be included in the recovery study, though it was not possible to find precursor and product ions for the two former compounds and dithianon was found to be lost during the extraction/clean-up procedure. Acibenzolar acid has though been reported in the literature to be possible to analyse by LC-MS/MS in negative ionization mode. However as it is an acid clean-up using PSA should probably be avoided. Chloropicrin has also by others been reported to be difficult to analyse by both GC-MS/MS and LC-MS/MS. Dithianon is sensitive to heat and may therefore be degraded during the extraction because of the increased temperatures resulting from the addition of  $\text{MgSO}_4$ .

### Overall validation on all 3 cereal types.

Of the 34 compounds included in the validation study (Table 1) 31 compounds were successfully validated on all three cereal types analysed by GC-MS/MS (10 pesticides), LC-MS/MS (14 pesticides) or both (7 pesticides), see **Appendix 2**.

For the accepted validation parameters the relative repeatability ( $\text{RSD}_r$ ) varied between 3-18 % with an average on 8%. The internal reproducibility ( $\text{RSD}_R$ ) varied between 5-23% with an average on 12%. Recoveries was in the range of 50-123% at all three concentration levels with an average on 88%. The combined LOQs were in the range of 0.01-0.02 mg/kg. Recoveries down to 50% was accepted if similar for all three spike levels and  $\text{RSD}_r$  and  $\text{RSD}_R$  were relatively low.

Asulam, foramsulfuron and TFNA were not possible to validate neither on GC-MS/MS or LC-MS/MS. For asulam the recovery were low (<50%) and at the two lowest spike levels were  $\text{RSD}_r$  and  $\text{RSD}_R$  >20%. For foramsulfuron low recoveries were found, though the recovery was equally low for all three spike levels and  $\text{RSD}_r$  and  $\text{RSD}_R$  were <20%. Foramsulfuron belong to the group of sulfonylureas and most of these are unstable at low pH and may not be protected sufficiently at the pH of the extract. Further if extract are stored for several days during analysis and perhaps reanalysis the compound may be degraded. Foramsulfuron is also a weak acid why some may be bound to PSA during d-SPE. TFNA resulted in an insufficient response by LC-MS/MS to produce data useful for validation and were not GC-MS/MS amenable. TFNA which is also an acid is most likely bound to PSA during the d-SPE. TFNA may also on some instruments give a better signal in ESI negative mode (for more information refer to : [http://www.crl-pesticides.eu/userfiles/file/EurlSRM/EurlSRM\\_meth\\_FlonicamidMetabolites.pdf](http://www.crl-pesticides.eu/userfiles/file/EurlSRM/EurlSRM_meth_FlonicamidMetabolites.pdf))

### **Validation on individual cereal type.**

More or less similar results were obtained if calculating the validation parameters for each of the cereal types individually.

Though an LOQ of 0.1 could be obtained for metobromuron if only including the data obtained for one cereal type whereas it was not possible to validate in the overall validation due to  $RSD_R > 20\%$ .

The validation results obtained for the individual cereals types are presented in Appendix 3 (wheat), 4 (barley) and 5 (rice).

## **8. Conclusions**

In conclusion 31 pesticides were successfully validated on wheat, barley and rice using the QuEChERS method and GC-MS/MS or/and LC-MS/MS. The LOQ obtained were 0.01 mg/kg except for three pesticides for which an LOQ of 0.02 mg/kg was obtained.

## **9. References**

- 1** EN 15662:2008. Foods of plant origin - Determination of pesticide residues using GC-MS and/or LC-MS/MS following acetonitrile extraction/partitioning and clean-up by dispersive SPE - QuEChERS-method
- 2** ISO 5725-2:1994. Accuracy (trueness and precision) of measurement methods and results – Part2. Basic method for the determination of repeatability and reproducibility of standard measurement method. First edition. December 1994.
- 3** Method Validation and Quality Control Procedures for Pesticide Residue Analysis in Food and Feed, Document No SANCO/12495/2011, 01/01/2012, European Commission, Brussels, 2012.

**Appendix 1a. MRM transitions GC-MS/MS.**

<b>GC-MS/MS</b>	<b>Retention time</b>	<b>Precursor ion-1</b>	<b>Product ion-1</b>	<b>CE</b>	<b>Precursor ion-2</b>	<b>Product ion-2</b>	<b>CE</b>
Acibenzolar-S-Methyl	13.1	182	181	5	182	153	22
Anthraquinone	14.1	180	152	15	208	180	10
Beflubutamid	15.5	221	193	10	355	176	10
Bifenazate	21.5	300	258	10	300	196	15
Biphenyl	8.0	153	152	15	154	153	15
Butralin	14.7	266	220	15	266	190	15
Carfentrazone-Ethyl	19.4	411	340	10	411	330	10
Dichlobenil	7.8	173	138	15	171	136	15
Dimethenamid	12.7	232	154	10	230	154	10
Ethalfuralin	10.0	276	202	10	316	276	10
Ethofumesate	13.7	207	161	10	286	207	12
Etoxazole	21.8	300	270	22	302	274	15
Flumioxazin	28.2	354	326	10	354	312	10
Isocarbophos	14.4	230	212	10	136	108	15
Propachlor	9.7	196	120	10	176	120	10
Pyraflufen-Ethyl	19.9	412	349	15	349	307	15
Tetramethrin	21.5	164	135	10	164	107	17

**Appendix 1b. MRM transitions for LC-MS/MS.**

LC-MS/MS	Retenti on time	Precursor ion-1	Product ion-1	CV	CE	Precursor ion-2	Product ion-2	CV	CE
2-hydroxypropoxy-carbazone 4-trifluoromethyl-nicotinic acid (TNFA)	8.1	415.1	174	40	12	415.1	116	40	10
Asulam	6.7	192.0	78.8	40	30	192.0	97.9	40	26
AvermectinB1a	7.3	230.9	156	40	20	230.9	92	40	20
Butralin	25.5	890.6	567.4	40	10	890.6	305.1	40	30
Carfentrazone-Ethyl	24.0	296.0	240	30	12	296.0	222	30	20
Chlorotoluron	19.3	412.2	346	40	20	412.2	366	40	17
Dimethenamid	14.3	213.0	72	20	20	213.0	140	20	30
Ethofumesate	16.4	276.0	244	40	10	276.0	168	40	20
Ethoxysulfuron	16.3	287.0	258.9	30	12	287.0	162.2	30	16
Fenpyroximate (E) (Z)	17.1	399.1	261	40	16	399.1	218	40	28
Foramsulfuron	24.5	422.2	366	50	10	422.2	135	50	30
Forchlorfenuron	13.4	453.2	182	20	23	453.2	254	20	20
Imazosulfuron	15.6	248.0	129	40	18	248.0	155	40	14
Meptyldinocap (ESI-)	16.4	413.0	156	30	20	413.0	153	30	10
Metobromuron	24.8	295.2	192.8	20	40	295.2	163	20	40
Nicosulfuron	14.5	258.9	170	40	30	258.9	148	40	20
Propachlor	12.0	411.1	182	30	20	411.1	212.9	30	20
Propoxycarbazone-sodium salt	14.5	211.9	170	30	10	211.9	94	30	27
Pyraflufen-Ethyl	9.9	421.0	180	20	10	421.0	137.9	20	27
Rotenone	19.7	413.0	338.9	30	17	413.0	289	30	30
Spirotetramat	18.9	395.1	213	20	26	395.1	192	20	26
Tetramethrin	17.8	374.2	330.2	20	10	374.2	302.2	20	20
Trinexapac-ethyl	22.3	332.2	164	40	23	332.2	135	40	17
	14.8	253.0	207	20	12	253.0	69	20	12

**Appendix 2. Recoveries, repeatability (RSD<sub>r</sub>), internal reproducibility (RSD<sub>R</sub>) and Limit of Quantification (LOQ) for pesticides validated on three cereal commodities, wheat, barley and rice.**

*Numbers in italic is outside 70-120% recovery or above 20% RSD*

	Wheat, barley and rice - QuEChERS	Spike level 0.01 mg/kg			Spike level 0.02 mg/kg			Spike level 0.1 mg/kg			LOQ
		Recovery %	RSD <sub>r</sub> , %	RSD <sub>R</sub> , %	Recovery %	RSD <sub>r</sub> , %	RSD <sub>R</sub> , %	Recovery %	RSD <sub>r</sub> , %	RSD <sub>R</sub> , %	
LC	2-hydroxypropoxycarbazone	71	13	16	68	9	8	67	6	8	0.01
GC	Acibenzolar-S-methyl	91	7	23	94	4	18	95	3	14	0.01
GC	Anthraquinone	88	7	11	96	5	13	93	6	12	0.01
LC	Avermectin B1a	99	16	15	102	16	16	96	10	11	0.01
GC	Beflubutamid	90	7	8	99	4	9	100	4	9	0.01
GC	Bifenazate	75	7	18	73	6	14	72	4	10	0.01
GC	Biphenyl		a, b, f		123	9	11	99	6	13	0.02
LC	Butralin	78	18	16	81	11	12	93	6	10	0.01
GC	Butralin	90	10	22	88	5	16	87	6	15	0.01
LC	Carfentrazone-ethyl	75	15	23	88	15	21	96	7	10	0.01
GC	Carfentrazone-ethyl	83	5	14	90	4	12	95	4	8	0.01
LC	Chlorotoluron	89	8	9	88	5	11	88	4	10	0.01
GC	Dichlobenil	98	9	16	102	6	10	102	5	13	0.01
LC	Dimethenamid	88	6	11	87	3	9	90	4	5	0.01
GC	Dimethenamid	94	5	12	99	6	12	97	5	11	0.01
GC	Ethalfuralin	95	12	12	96	13	15	97	9	12	0.01
LC	Ethofumesate	79	15	15	91	5	6	94	5	6	0.01
GC	Ethofumesate	94	8	7	103	7	7	103	5	6	0.01
LC	Ethoxysulfuron	75	9	16	75	6	13	80	5	8	0.01
GC	Etoxazole	92	8	14	94	7	12	95	5	11	0.01
LC	Fenpyroximate	86	6	11	90	4	10	95	5	8	0.01
GC	Flumioxazin	93	9	21	98	7	15	96	6	10	0.01
LC	Forchlorfenuron	77	13	23	75	7	13	80	5	7	0.01

Wheat, barley and rice - QuEChERS		Spike level 0.01 mg/kg			Spike level 0.02 mg/kg			Spike level 0.1 mg/kg			
		Recovery %	RSD <sub>r</sub> , %	RSD <sub>R</sub> , %	Recovery %	RSD <sub>r</sub> , %	RSD <sub>R</sub> , %	Recovery %	RSD <sub>r</sub> , %	RSD <sub>R</sub> , %	LOQ
LC	Imazosulfuron	68	12	12	71	7	8	74	6	7	0.01
GC	Isocarbophos	96	9	22	104	8	18	101	4	15	0.01
LC	Meptyldinocap		a, b, j		95	17	16	102	14	13	0.02
LC	Metobromuron	83	14	22	82	13	18	85	5	9	0.01
LC	Nicosulfuron	58	10	11	55	6	9	55	7	8	0.01
LC	Propachlor	86	8	11	90	4	9	93	4	10	0.01
GC	Propachlor		b		98	10	10	98	7	11	0.02
LC	Propoxycarbazone sodium salt	81	5	8	79	7	7	80	5	6	0.01
LC	Pyraflufen-ethyl	82	7	8	86	4	11	91	4	11	0.01
GC	Pyraflufen-ethyl	88	6	9	96	5	7	99	4	7	0.01
LC	Rotenone	85	11	13	90	6	12	93	9	13	0.01
LC	Spirotetramat	79	8	12	85	4	11	91	4	8	0.01
LC	Tetramethrin	86	8	8	91	5	10	98	5	10	0.01
GC	Tetramethrin	95	13	20	102	8	16	95	6	12	0.01
LC	Trineexpac-ethyl	58	13	17	52	14	19	50	6	12	0.01
<b>Pesticides/metabolites not possible to validate</b>											
LC	Asulam		a, b, e								
LC	Foramsulfuron		e								
LC	TFNA		h								

a) RSD<sub>r</sub> > 20%; b) RSD<sub>R</sub> > 20%; c) Not GC-MS/MS amenable; d) Not LC-MS/MS amenable; e) Recovery <50%; f) Recovery >50%; h) not multimetod amenable; j) wheat results not included because not possible to validate for this matrix.

**Appendix 3. Recoveries, repeatability (RSD<sub>r</sub>) and Limit of Quantification (LOQs) for pesticides validated on wheat.**

	Wheat - QuEChERS	Spike level, mg/kg 0.01		Spike level, mg/kg 0.02		Spike level, mg/kg 0.1		LOQ
		Recovery %	RSDr %	Recovery %	RSDr %	Recovery %	RSDr %	
LC	2-hydroxypropoxycarbazone	63	16	68	15	63	8	0.01
GC	Acibenzolar-S-methyl	104	2	105	5	105	3	0.01
GC	Anthraquinone	96	5	106	4	101	9	0.01
LC	Avermectin B1a	96	15	111	12	102	9	0.01
GC	Beflubutamid	95	8	107	6	104	4	0.01
GC	Bifenazate	82	8	81	5	80	3	0.01
GC	Biphenyl	91	23	113	7	112	7	0.02
GC	Butralin	95	7	93	5	93	5	0.01
LC	Butralin	78	9	82	13	86	7	0.01
LC	Carfentrazone ethyl	65	18	71	18	87	8	0.01
GC	Carfentrazone-ethyl	73	4	81	3	94	6	0.01
LC	Chlorotoluron	90	9	89	2	84	4	0.01
GC	Dichlobenil	101	8	111	4	116	6	0.01
GC	Dimethenamid	106	5	110	5	107	6	0.01
LC	Dimethenamid	89	6	90	2	91	5	0.01
GC	Ethalfuralin	89	16	89	17	98	14	0.01
GC	Ethofumesate	92	5	106	6	105	7	0.01
LC	Ethofumesate	83	15	94	5	97	5	0.01
LC	Ethoxysulfuron	85	7	85	5	86	6	0.01
GC	Etoxazole	101	8	102	7	101	7	0.01
LC	Fenpyroximate	90	5	90	4	92	6	0.01
GC	Flumioxazin	107	3	111	4	105	8	0.01
LC	Foramsulfuron	52	16	49	13	47	9	0.01
LC	Forchlorfenuron	82	10	80	4	84	6	0.01
LC	Imazosulfuron	68	11	73	11	74	7	0.01



Wheat - QuEChERS		Spike level, mg/kg 0.01		Spike level, mg/kg 0.02		Spike level, mg/kg 0.1		
		Recovery %	RSDr %	Recovery %	RSDr %	Recovery %	RSDr %	LOQ
GC	Isocarbophos	101	9	104	6	103	4	0.01
GC	Metobromuron	a, e, g		a, e, g		116	7	0.1
LC	Metobromuron	78	12	85	6	91	5	0.01
LC	Nicosulfuron	60	11	58	2	54	9	0.01
GC	Propachlor	87	11	98	5	108	6	0.01
LC	Propachlor	88	7	89	4	87	4	0.01
LC	Propoxycarbazone sodium salt	87	4	79	6	76	7	0.01
LC	Pyraflufen ethyl	81	9	79	3	81	5	0.01
GC	Pyraflufen-ethyl	86	6	97	3	103	5	0.01
LC	Rotenone	a		a		82	16	0.1
LC	Spirotetramat	73	8	79	3	86	5	0.01
GC	Tetramethrin	105	9	111	4	105	7	0.01
LC	Tetramethrin	87	10	88	5	91	5	0.01
<b>Pesticides/metabolites not possible to validate</b>								
LC	Asulam	a, e						
LC	Bifenazate	e						
LC	Meptyldinocap	i						
LC	TFNA	h						
LC	Trineexpac-ethyl	e						

a) RSDr > 20%; e) Recovery <50%; g) To low sensitivity; h) not multimetod amenable ; i) interfering matrix enables quantification.

#### Appendix 4. Recoveries, repeatability (RSD<sub>r</sub>) and Limit of Quantification (LOQs) for pesticides validated on barley.

Barley - QuEChERS		Spike level, mg/kg 0.01		Spike level, mg/kg 0.02		Spike level, mg/kg 0.1		
		Recovery %	RSD <sub>r</sub> %	Recovery %	RSD <sub>r</sub> %	Recovery %	RSD <sub>r</sub> %	LOQ
LC	2-hydroxypropoxycarbazone	78	15	66	4	67	5	0.01
GC	Acibenzolar-S-methyl	67	8	75	6	80	3	0.01
GC	Antraquinone	79	9	83	7	82	5	0.01
LC	Asulam	70	19	62	12	55	10	0.01
LC	Avermectin B1a	a		62	16	93	5	0.02
GC	Beflubutamid	85	6	92	4	91	5	0.01
GC	Bifenazate	61	8	62	7	68	5	0.01
GC	Biphenyl	212	22	132	9	88	6	0.02
GC	Butralin	69	17	73	7	73	8	0.01
LC	Butralin	a		75	15	91	6	0.02
LC	Carfentrazone ethyl	a		96	20	99	6	0.02
GC	Carfentrazone-ethyl	81	5	89	5	88	5	0.01
LC	Chlorotoluron	8	9	81	4	82	4	0.01
GC	Dichlobenil	83	8	94	10	92	6	0.01
GC	Dimethenamid	85	6	89	11	87	6	0.01
LC	Dimethenamid	79	8	79	5	87	3	0.01
GC	Ethalfuralin	97	10	96	15	89	7	0.01
GC	Ethofumesate	93	10	98	7	98	5	0.01
LC	Ethofumesate	77	17	88	5	91	6	0.01
LC	Ethoxysulfuron	75	10	74	8	81	4	0.01
GC	Etoxazole	80	8	83	7	85	5	0.01
LC	Fenpyroximate	77	8	82	2	91	4	0.01
GC	Flumioxazin	72	14	83	9	88	5	0.01
LC	Forchlorfenuron	60	16	66	11	77	5	0.01
LC	Imazosulfuron	71	8	73	3	76	6	0.01

Barley - QuEChERS		Spike level, mg/kg 0.01		Spike level, mg/kg 0.02		Spike level, mg/kg 0.1		LOQ
		Recovery %	RSDr %	Recovery %	RSDr %	Recovery %	RSDr %	
GC	Isocarbophos	75	7	86	10	86	6	0.01
LC	Meptyldinocap	119	24	96	18	102	18	0.01
GC	Metobromuron	a, e, g		a, e, g		71	7	0.1
LC	Metobromuron	72	19	71	14	79	3	0.01
LC	Nicosulfuron	60	12	55	6	58	3	0.01
GC	Propachlor	85	17	94	15	93	8	0.01
LC	Propachlor	79	12	84	5	90	3	0.01
LC	Propoxycarbazone sodium salt	78	8	77	8	81	3	0.01
LC	Pyraflufen ethyl	79	9	83	5	93	2	0.01
GC	Pyraflufen-ethyl	83	5	90	5	91	5	0.01
LC	Rotenone	81	14	84	7	96	5	0.01
LC	Spirotetramat	75	9	81	7	90	4	0.01
GC	Tetramethrin	78	17	86	11	85	6	0.01
LC	Tetramethrin	83	8	83	3	95	6	0.01
LC	Trineexpac-ethyl	a		55	6	52	5	0.02
<b>Pesticides/metabolites not possible to validate</b>								
LC	Asulam	a, e						
LC	Foramsulfuron	e						
LC	TFNA	h						

a) RSDr > 20%; e) Recovery < 50%; g) To low sensitivity; h) not multimethod amenable.

**Appendix 5. Recoveries, repeatability (RSD<sub>r</sub>) and Limit of Quantification (LOQs) for pesticides validated on rice.**

	Rice - QuEChERS	Spike level, mg/kg 0.01		Spike level, mg/kg 0.02		Spike level, mg/kg 0.1		LOQ
		Recovery %	RSDr %	Recovery %	RSDr %	Recovery %	RSDr %	
LC	2-hydroxypropoxycarbazone	71	8	69	7	72	4	0.01
GC	Acibenzolar-S-methyl	101	9	101	2	99	4	0.01
GC	Anthraquinone	90	7	98	3	95	4	0.01
LC	Avermectin B1a	101	16	102	13	94	15	0.01
GC	Beflubutamid	89	5	98	2	105	2	0.01
GC	Bifenazate	83	5	76	6	69	2	0.01
LC	Bifenazate	63	10	66	9	61	6	0.01
GC	Biphenyl	148	17	124	10	99	3	0.02
GC	Butralin	104	7	98	5	95	6	0.01
LC	Butralin	78	24	87	8	102	6	0.01
LC	Carfentrazone ethyl	84	13	96	8	100	7	0.01
GC	Carfentrazone-ethyl	95	4	101	4	102	3	0.01
LC	Chlorotoluron	93	6	96	4	97	4	0.01
GC	Dichlobenil	109	10	102	4	99	4	0.01
GC	Dimethenamid	93	6	97	3	98	2	0.01
LC	Dimethenamid	95	3	93	2	92	3	0.01
GC	Ethalfuralin	99	9	104	7	105	5	0.01
GC	Ethofumesate	96	6	104	8	105	3	0.01
LC	Ethofumesate	78	14	91	6	94	2	0.01
LC	Ethoxysulfuron	65	9	67	6	75	4	0.01
GC	Etoxazole	95	9	97	6	100	3	0.01
LC	Fenpyroximate	92	5	98	4	103	4	0.01
GC	Flumioxazin	99	9	100	8	96	4	0.01
LC	Forchlorfenuron	88	13	82	6	78	5	0.01
LC	Imazosulfuron	65	17	66	5	70	6	0.01

Rice - QuEChERS		Spike level, mg/kg 0.01		Spike level, mg/kg 0.02		Spike level, mg/kg 0.1		LOQ
		Recovery %	RSDr %	Recovery %	RSDr %	Recovery %	RSDr %	
GC	Isocarbophos	113	10	121	7	115	3	0.01
LC	Meptyldinocap		a	95	19	103	9	0.02
GC	Metobromuron	a, e, g		a, e, g		98	14	0.1
LC	Metobromuron	100	11	87	7	86	6	0.01
LC	Nicosulfuron	53	1	50	9	53	7	0.01
GC	Propachlor	119	20	102	9	94	6	0.01
LC	Propachlor	92	5	98	4	104	4	0.01
LC	Propoxycarbazone sodium salt	78	3	81	7	82	5	0.01
LC	Pyraflufen ethyl	86	4	96	5	99	4	0.01
GC	Pyraflufen-ethyl	95	6	101	5	102	3	0.01
LC	Rotenone	90	9	97	5	102	3	0.01
LC	Spirotetramat	87	7	94	4	98	4	0.01
GC	Tetramethrin	101	15	108	9	96	5	0.01
LC	Tetramethrin	89	5	99	6	107	3	0.01
LC	Trineexpac-ethyl	66	5	54	9	55	3	0.01
<b>Pesticides/metabolites not possible to validate</b>								
LC	Asulam	a, e						
LC	Foramsulfuron	e						
LC	TFNA	h						

a) RSDr > 20%; b) RSDR > 20%; c) Not GC-MS/MS amenable; d) Not LC-MS/MS amenable; e) Recovery <50%; f) Recovery >50%; g) To low sensitivity; h) not multimetod amenable; i) interfering matrix enables quantification

## Appendix 6: Principles of the QuEChERS method for cereal extraction

# QuEChERS for cereals (FP417)

Weigh 5 g ( $\pm 0.05$  g) of flour into a 50 ml single use centrifuge tube (red cap).  
Add internal standard and/or spike standard (maximum 25  $\mu$ l)

Add a ceramic homogenizer and 10 g of cold water and shake briefly

Add 10 ml acetonitrile and shake vigorously by hand for 1 min. (1. extraction)

Add the prepared mixture of 4 g  $\text{MgSO}_4$ , 1 g NaCl, 1 g  $\text{Na}_3$  citrate dihydrate and 0.5 g  $\text{Na}_2\text{H}$  citrate sesquihydrate. Shake for a few seconds after each addition to prevent lumps.

Shake vigorously for 1 min. (2. Extraction with phase separation)

Centrifuge for 10 min at 4500 rpm

Transfer at least 8 ml of the extract to a 15 ml single use centrifuge tube and store in the freezer ( $-80^\circ\text{C}$  for 1 hour or over night). When the extract are almost thawed (i.e. About  $-40^\circ\text{C}$ ) centrifugate (should be cold  $5^\circ\text{C}$ ) for 5 min. at 4500 rpm.

Transfer 6 ml of the cold extract to a 15 ml single use centrifuge tube containing 150 mg PSA and 900 mg  $\text{MgSO}_4$ . Close the tube and shake vigorously for 30 seconds.

Centrifuge for 5 min. at 4500 rpm

Transfer 4 ml of the extract to a 15 ml single use centrifuge tube. Add 40  $\mu$ l of 5% formic acid solution in acetonitrile (10  $\mu$ l/ml extract). Dilute the extract 1:1 with acetonitrile

Transfer the final extract into auto sampler vials and analyse by GC and LC.